


Fluoromount-G[®] Anti-Fade Mounting Medium

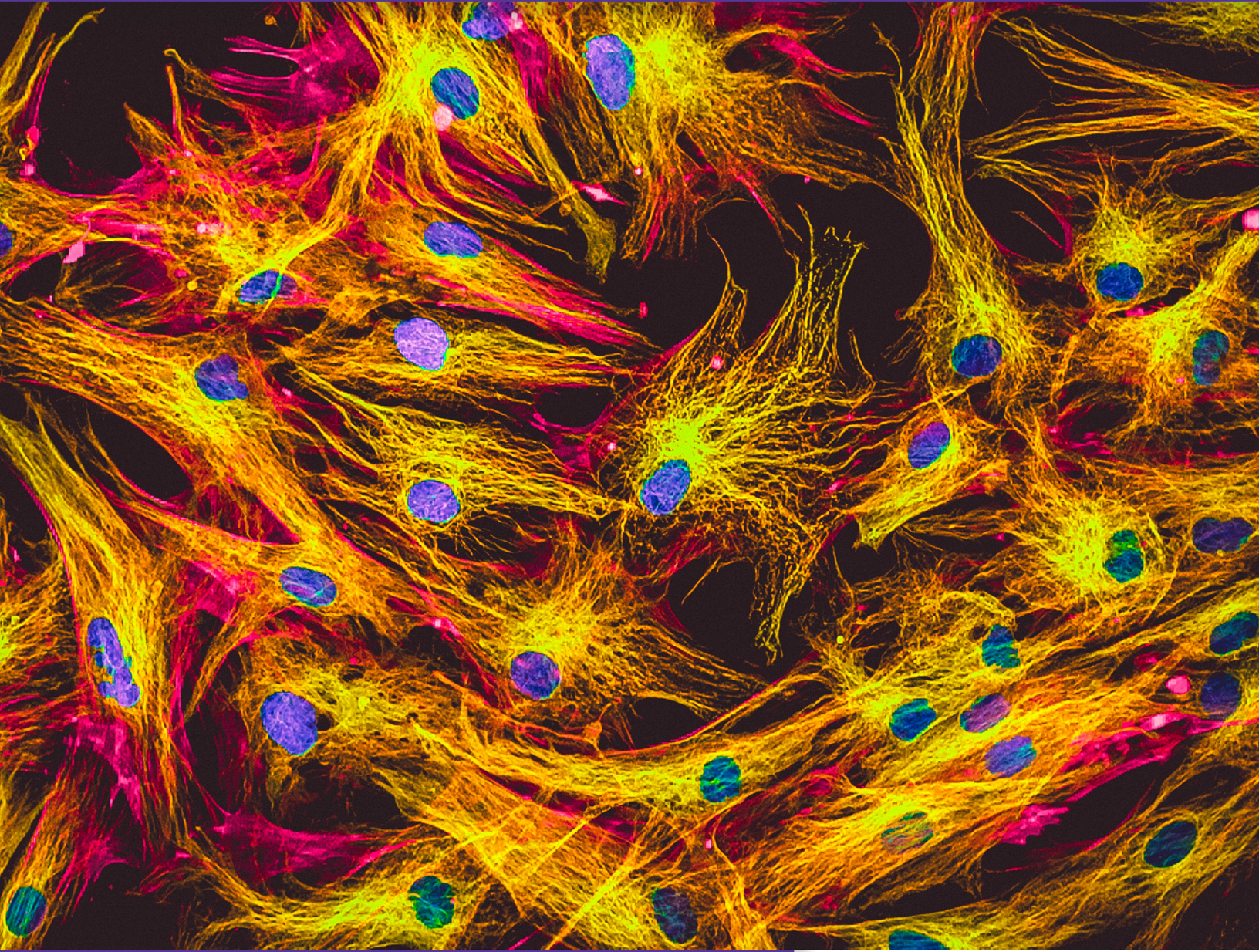

Quality Antibodies for Quality Research

SouthernBiotech



www.SouthernBiotech.com info@SouthernBiotech.com




Fluoromount-G® Anti-Fade Mounting Medium

Mounting Medium	Cat. No.	Qty
Fluoromount-G® Anti-Fade	0100-35	20 mL

Fluoromount-G® Anti-Fade is an aqueous-based compound for mounting slides used in microscopy applications in which the preservation of fluorescence intensity is critical. This mounting medium exhibits minimal background in the blue spectrum and provides a semi-permanent seal for storage of slide preparations.

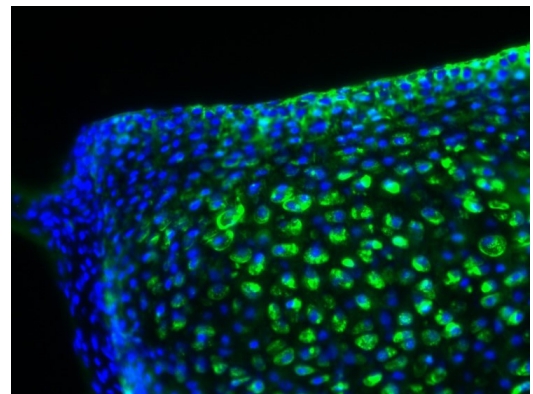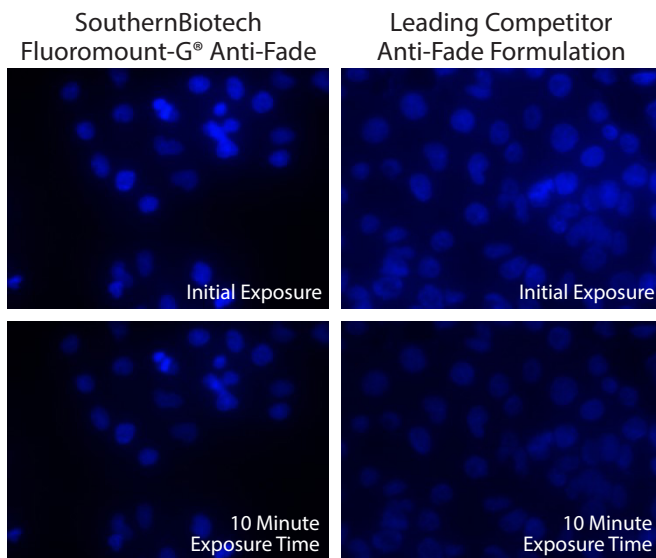

Frozen newborn mouse cartilage section was blocked with Normal Rabbit Serum (SB Cat. No. 0040-01) and stained with Mouse Anti-Collagen Type II-UNLB followed by Rabbit F(ab')₂ Anti-Mouse IgG(H+L), Human ads-BIOT (SB Cat. No. 6125-08), Streptavidin-FITC (SB Cat. No. 7100-02), DAPI, and mounted with Fluoromount-G® Anti-Fade (SB Cat. No. 0100-35).

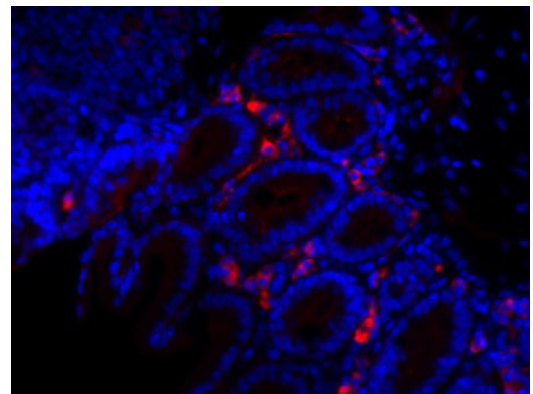

Paraffin embedded mouse small intestine tissue section was blocked with Normal Rabbit Serum (SB Cat. No. 0040-01) and stained with Goat Anti-Mouse IgA-UNLB (SB Cat. No. 1040-01) followed by Rabbit Anti-Goat IgG(H+L), Human SP ads-BIOT (SB Cat. No. 6164-08), Streptavidin-CY3.5 (SB Cat. No. 7100-24), DAPI, and mounted with Fluoromount-G® Anti-Fade (SB Cat. No. 0100-35).

Corporate Offices

160 Oxmoor Blvd
Birmingham, AL 35209
USA

Mailing Address

PO Box 26221
Birmingham, AL 35260
USA

Tel 205.945.1774

Fax 205.945.8768

Email info@SouthernBiotech.com

Technical Support technical@SouthernBiotech.com

Website www.SouthernBiotech.com

SouthernBiotech is ISO 9001:2015 Certified

SouthernBiotech

www.SouthernBiotech.com info@SouthernBiotech.com



Fluoromount-G® is a registered trademark of Southern Biotechnology Associates, Inc.
Cy® is a registered trademark of GE Healthcare
MM-FAF-00