

Buffers and Stabilizers



PROVIDING CONSCIOUS IMMUNOASSAY SOLUTIONS

kementec


Buffers and Stabilizers

Kementec Solutions' ECO-TEK buffers are designed for use in colorimetric as well as in chemiluminescence based assays such as ELISA, Western Blotting and IHC. They are all ready-to-use and robust diluents offering superior stabilization of proteins conjugated to Horseradish Peroxidase (HRP) and Alkaline Phosphatase (AP). These stabilizers prevent loss of enzyme activity and maintain the proteins' molecular conformation. Furthermore all our ECO-TEK buffers are Bovine Serum Albumin (BSA) free and several are animal free.

WellChampion®

Plate Blocking & Protection Buffer for pre-coated plates

WellChampion preserves and protects biological activity of proteins and other biomolecules coated and dried on microwell surfaces. The protection of the coated and dried surface prevents degradation, denaturation and leaching. As WellChampion also efficiently blocks any free binding sites left on the microwell surface after coating, separate blocking of the surface is not necessary. WellChampion dramatically speeds up your plate coating and stabilization process in most cases at the same time resulting in improved assay sensitivity.

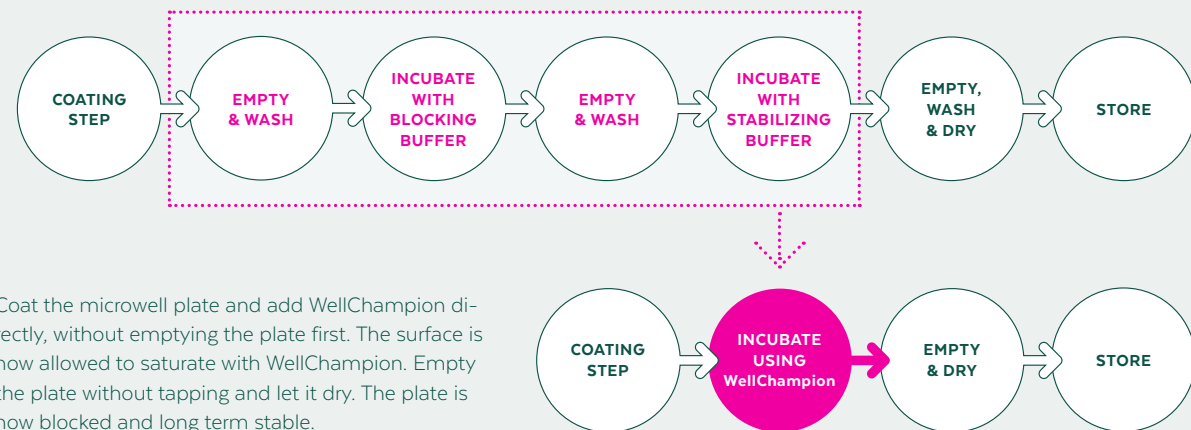

Coat the microwell plate and add WellChampion directly, without emptying the plate first. The surface is now allowed to saturate with WellChampion. Empty the plate without tapping and let it dry. The plate is now blocked and long term stable.

Choose the right buffer

The following table will guide you to the optimal choice of our different solutions.
D: Diluent, S: Stabilizer, B: Blocker

Animal-free Buffers	Standards & Controls	Samples	HRP-conjugates	AP-conjugates	Antibodies	Surfaces
HRP-StabilPLUS *			D/S			
AP-StabilPLUS				D/S		
UNI-StabilPLUS			D/S	D/S		
Protein-StabilPLUS	D/S				D/S	
Antibody Enhancer *			D/S	D/S	D/S	
Sample X-tra	D	D				
Effect Diluent		D	D	D	D	
Synthetic Blocker						B
WellChampion® *						B/S

* Not Animal-free

Synthetic Blocking Buffers



Blocking Buffers for ELISA and Blotting

Purely synthetic, protein free blocking buffers for ELISA and Blotting applications with high blocking capacity: Improving signal-to-noise ratio values and eliminating the risk of false positive cross reactions in the assay.

HRP-StabilPLUS

HRP Conjugate Diluent & Stabilizer

HRP-StabilPLUS is our most successful and widely used conjugate buffer for HRP labeled proteins due to its unsurpassed long-lasting stabilizing effect in both high and low protein concentrations.

AP-StabilPLUS



AP Conjugate Diluent & Stabilizer

AP-StabilPLUS is a highly recommended alternative to conventional Alkaline Phosphatase conjugate stabilizers. AP-StabilPLUS provides a long-lasting stabilizing effect.

UNI-StabilPLUS



AP & HRP Conjugate Diluent & Stabilizer

UNI-StabilPLUS is a universal conjugate diluent designed for the stabilization of both AP- & HRP-labelled proteins. Utilizing the same buffer for all conjugates simplifies daily work processes.

Antibody Enhancer

HRP- & AP- Antibody Diluent, Stabilizer & Amplifier

Antibody Enhancer is an antibody efficient conjugate stabilizer as well as signal amplifier. The buffer enhances the binding capacity of the diluted antibodies and provides the highest assay signal-to-noise ratio among our stabilizers.

Sample X-tra



Sample dilution buffer for Immunoassays

The dilution of sample matrices & assay standards with Sample X-tra minimizes nonspecific interactions between the matrix proteins and the surface of the plate, improving signal-to-noise ratio.

Effect Diluent buffers, Animal-free



Sample dilution buffers for Immunoassays

The 3 Effect Diluent buffers (Low, Medium and High) eliminate interference from HAMAs (human-a-mouse Abs), Heterophilic Abs and RF (rheumatoid factors) by blocking cross reactions without affecting the correct binding of the assay antibodies.

Protein-StabilPLUS



Protein Stabilizer

Protein-StabilPLUS gives superior stabilization of ready-to-use pre-diluted standards and controls used in immunoassays.



Animal-free



Our immunoassay solutions are developed in order to create a healthy and safe work environment and preserve natural resources, while helping our customers fulfill significant regulatory requirements.

We call it ECO-TEK®

Profile

Kementec is a well-established, innovative, and service oriented company with more than 35 years of experience as a manufacturer and supplier of components to the industry and to researchers worldwide.

Being a market leader in developing solutions for immunoassays, our continuous goal is to deliver sustainable products without compromising on quality, sensitivity and accuracy.

Kementec values long lasting business relationships and the success of our customers. We consider ourselves more as a trusted partner than as a general supplier. A part of this vision is providing high quality colorimetric substrates and color-coded reagents, described as:

Colors of Success!

Kementec Solutions A/S
Taastrup, Denmark

Tel. +45 39 27 17 71
support@kementec.com

Kementec Solutions Inc.
Amherst, NH

Tel. (617) 336-7474
contact@kementec.com

kementec.com

kementec