



Biopanda Reagents is an ISO13485 certified company which designs, manufactures, and supplies in vitro test kits for clinical laboratories, veterinary practice, environmental protection, and food safety. Our company provides innovative, high quality and cost-effective test reagents that offer the sensitivity, specificity and reproducibility that customers demand.

Our range of veterinary products includes rapid test kits for detection of pathogens in feline, canine and bovine samples. These diagnostic kits utilise lateral flow technology, detecting pathogen specific antibodies present in an animal's bodily fluids or waste.

PRODUCT FOCUS:

Canine CRP Rapid Test, Product Number: RAPG-CCRP-001

Product Description

The Biopanda Canine CRP Rapid Test qualitatively detects Canine C-reactive protein (cCRP) from a dog's serum or plasma specimen.

Test Attributes


- Results available in only 5 minutes
- Easy to use
- Results can be read visually
- No need for an analyser
- Cost effective method for assisting to diagnose inflammation

Test Procedure

1. Collect a dog's fresh whole blood sample and separate serum or plasma from blood as soon as possible to avoid haemolysis.
2. Open the cassette from the foil pouch and position it on a flat, clean surface.
3. Holding the pipette vertically, slowly drip 2 drops of the serum or plasma sample into the specimen well and start the timer.
4. Read the results at 5 minutes.
After 10 minutes results are considered invalid.

Performance Characteristics

- Relative sensitivity: 92.86%
- Relative specificity: >99.9%

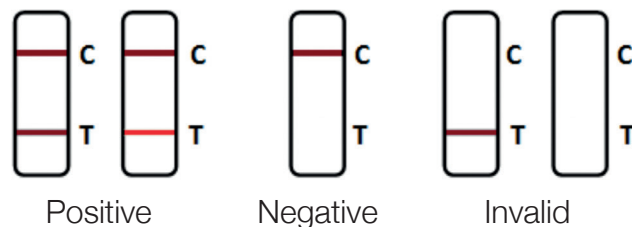

2BScientific is the exclusive UK distributor for Biopanda Reagents. For more information visit: www.2BScientific.com/Suppliers/Biopanda-Reagents

2BScientific
the life science company
with a difference

www.2BScientific.com

+44 (0)1869 238 033 sales@2BScientific.com